



Experts in Bottle Cleaning Machinery, Engineering and Fabrication Services

EXTERIOR VIAL WASHER/DRYER

The Exterior Vial Washer/Dryer will wash, rinse, and dry the exterior of filled and capped vials/ampoules. The entire exterior of each vial is exposed for cleaning as it passes through our patented “Zig-Zag” guide pattern. A low-pressure wash and rinse are used to keep the cap area dry. The option of multiple passes (1 to 10) including choice of rinse sequence makes the Exterior Vial Washer/Dryer the most versatile in the industry. Capable of reaching speeds of 150vpm - 500vpm (depending on vial size), this machine is the perfect addition to any pharmaceutical line.

Special Features

- Multiple rinses/dry passes allow for “bone dry” vials
- No Change Parts. Changeovers for the single pass washer are quick and easy, simply adjust the rails for new vials
- No cap wetting. The low-pressure wash gently cleans the exterior of vials, while leaving the caps spotless
- High Speed Cleaning. The 10-pass vial washer can reach speeds of up to 500vpm (depending on vial size)



E: sales@mcbrodyengineering.com

T: 815-744-8900



www.mcbrodyengineering.com

1251 S. Larkin Ave. Rockdale, IL 60436





Experts in Bottle Cleaning Machinery, Engineering and Fabrication Services

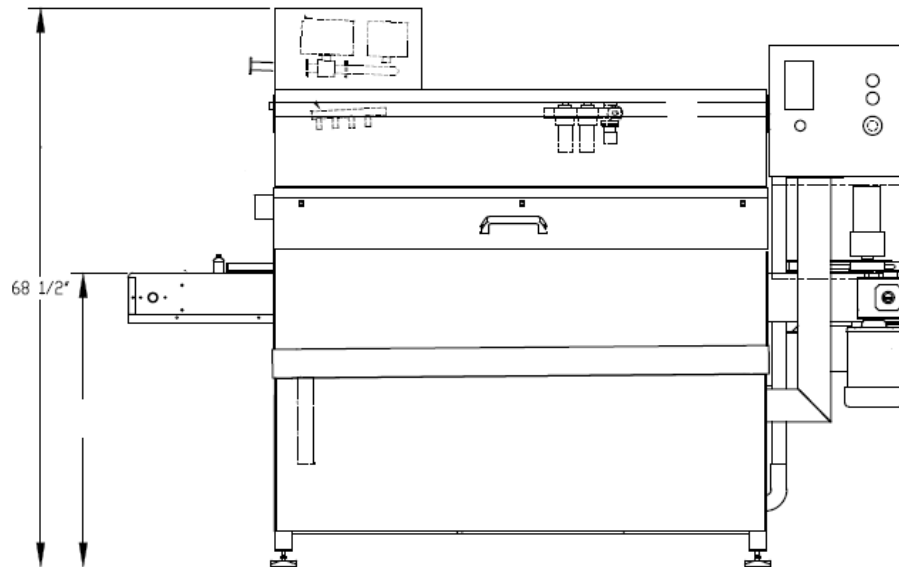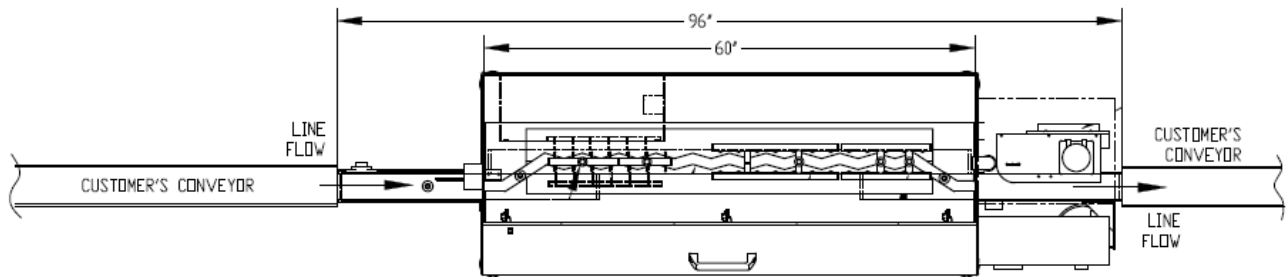

E: sales@mcbradyengineering.com

T: 815-744-8900

www.mcbradyengineering.com

1251 S. Larkin Ave. Rockdale, IL 60436

